

ANPEX

BY MELTON LAB

VAM-510

VETERINARY ANESTHESIA MACHINE

NEW MED



CE

Designed to be an excellent choice for you.

- ▶ Versatile breathing circuit
- ▶ Customizable configuration
- ▶ Qualified parts installed
- ▶ Interchangeable mounting methods
- ▶ User friendly design
- ▶ Safe oxygen flush


And more benefits need to be found ...

MELTON[™]
LABORATORY

Think more, do better.

ANPEX VAM-510

Veterinary Anesthesia Machine Series

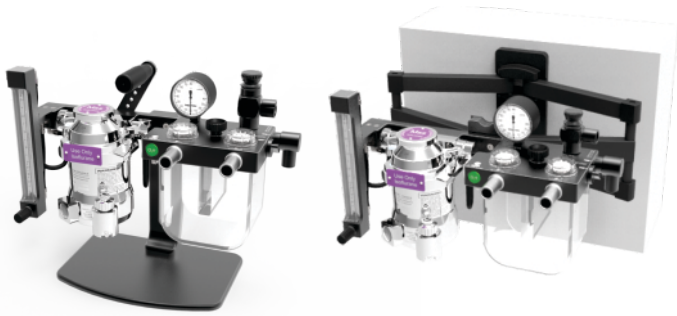

VAM-510*
Tabletop Type

VAM-510W
Wall Mount Type



VAM-510QS*
Quick Stand Mount

VAM-510S*
Stand Type

VAM-510ST*
Stand Type w/ Top Tray

* Indicates MRI-compatible types also available.

Benefits

Versatile Breathing Circuit

Opening or closing the rebreathing circuit by one turn of the knob, convenient to switch between different application scenarios without making changes to tubing system of the machine.

Safe Oxygen Flush

Oxygen flush valve with restricted flow for extra safety of the patients. Integrated with the breathing circuit for easy access and handle.

Interchangeable Mounting Methods

Tabletop design with easy carry handle, upgradable to stand-mount or wall-mount seamlessly by specific mounting kit.

Customizable Configuration

Flowmeter with different reading ranges, vaporizers of Cagemount and Selectatec type with screw-fill and key-fill option are available.

User Friendly Design

Delivery with a complete assembly, no tools needed for regular maintenance.

Qualified Parts Installed

All parts and materials used in ANPEX anesthesia machines are medical or food grade product from qualified global suppliers.

Specifications (VAM-510)

Dimension	452(L)x 208(W) x 360(H) mm
Weight (w/o vaporizer)	6.7 KG (14.8 lbs)
Flowmeter Reading	0.2 - 4LPM (Optional)
Manometer Reading	-40 ~ 80 cmH ₂ O
CO ₂ Absorber Capacity	1300 cc
APL Valve Adjustable Range	0 ~ 60 cmH ₂ O
Oxygen Flush Flowrate	10 ~ 15LPM
Oxygen Connector	D.I.S.S Male
Rated Working Pressure	50 PSIG (3.44 bar)

Accessories

We provide a wide range of high quality accessories for veterinarians along with our anesthesia machine. Contact us for a full catalogue and more detailed product specifications.

Anesthesia Accessories

- Disposable Breathing Circuits
- Non-rebreathing Circuits
- Laryngoscopes
- Endotracheal Tubes
- Animal Anesthesia Masks
- Breathing Bags



Gas Supply Components

- Medical Gas Regulator
- Pre-assembled Oxygen Hoses
- Hose Adaptors



Pollution Control

- F/Air Anesthesia Gas Filter Kit



Find more options by contacting with our sales representative or a distributor near you.

 sales@meltonlab.com


MELTON
LABORATORY


A PROUD MEMBER OF AALAS

MELTON LAB CO.,LTD
A SUBSIDIARY OF MELTON LAB LLC, TEMPLE CITY, CA 91780, USA

☎ +86.21.52271963

☎ Extension.200

📱 @MeltonLab

🌐 www.meltonlab.com